


Spring Tip Steel Guide Wire

Used to guide esophageal dilators for dilatation

Specially designed Flexible atraumatic tip
Spring base and metallic wire shaft are interlocked.
Joint of Spring and wire shaft is tapered to stop guide wire getting stuck in the dilator.
Smooth surface finish offers easy passage through dilator lumen during dilatation procedure.
Availability - Sterile, Single Use.

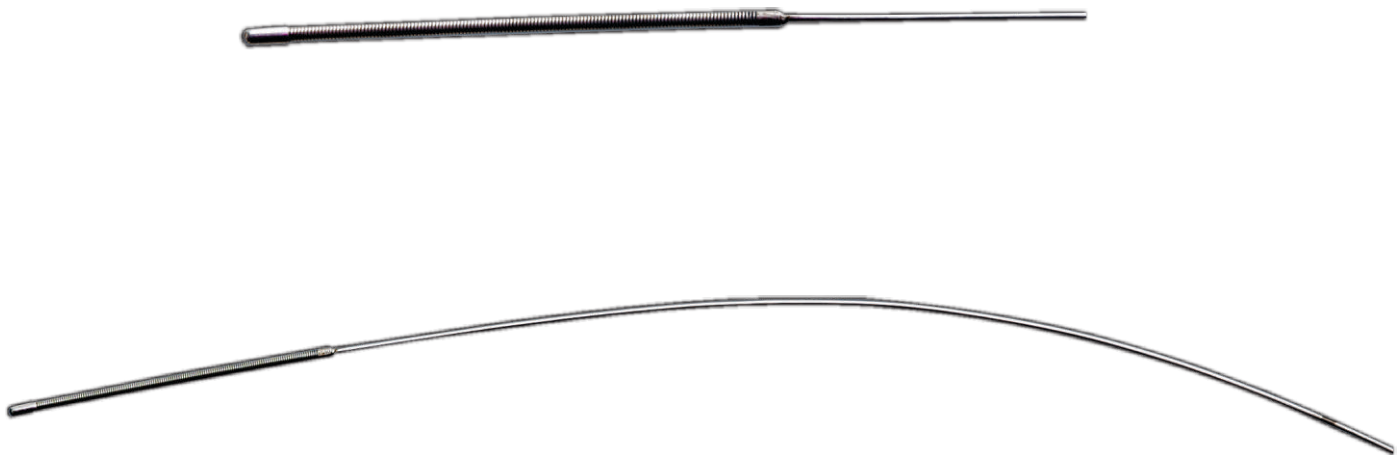

Ordering Information- Mention the product code x quantity required.

Reference	Product Name	Length	Tip Ø
SMGW	Spring Tip Steel Guide Wire	250 cm	2.00 mm
SMGW - M	Spring Tip Steel Guide Wire Marked	250 cm	2.00 mm

Tailormade device can be produced provided by pre approved drawing/specification